



IBA Lifesciences GmbH
Rudolf-Wissell-Str. 28
37079 Goettingen
Germany
Tel.: +49 (0) 551-5 06 72-0
E-mail: info@iba-lifesciences.com
www.iba-lifesciences.com

Data Sheet

pASG-IBA44

Cat. No.: 5-4044-001

Version: 3.0
Revision Date: 21.07.2021

Description	StarGate Acceptor Vector for bacterial expression. The expression cassette is under transcriptional control of the tetracycline promoter/operator. The expressed recombinant protein will be secreted into the periplasm.
Affinity tag	Strep-tag® II is fused to the N-terminus and 6xHis-tag is fused to the C-terminus of the recombinant protein.
Secretion	The ompA signal sequence directs the expressed protein into the periplasmic space and will be cleaved off during the translocation process.
Cloning Strategy	Cloning into StarGate Acceptor Vectors has to be done with the restriction enzyme Esp3I. There is no Multiple Cloning Site (MCS) available that can be used for the integration of the gene of interest instead (see manual).
Expression strain	Any <i>E. coli</i> strain. The <i>tet</i> -promoter works independently from the genetic background of <i>E. coli</i> .
Bacterial Expression	Expression is induced upon addition of 200 µg anhydrotetracycline per 1 liter <i>E. coli</i> shaking culture ($A_{550} = 0.5$).
Resistance	Ampicillin
Form	5 µg, dissolved in 20 µl TE buffer, pH 8.0: 10 mM Tris/HCl, 1 mM EDTA
Concentration	250 ng/µl
Stability	12 months after shipping
Storage	recommended: 2-8 °C for frequent usage, -20 °C for long-term storage
Shipping	room temperature
Hazards	Product is not classified as hazardous according to (EC) No 1272/2008 [CLP]. A Material Safety Data Sheet is provided.

For research use only


Trademark information

The owners of trademarks marked by "®" or "TM" are identified at <http://www.iba-lifesciences.com/patents.html>. Registered names, trademarks, etc. used in this document, even when not specifically marked as such, are not to be considered unprotected by law.

Important licensing information


This product is covered by intellectual property (IP) rights and on completion of the sale IBA Lifesciences grants respective Limited Use Label Licenses to purchaser. IP rights and Limited Use Label Licenses for said technology are further described and identified at <http://www.iba-lifesciences.com/patents.html> or upon inquiry at info@iba-lifesciences.com or at IBA Lifesciences GmbH, Rudolf-Wissell-Str. 28, 37079 Goettingen, Germany. By use of this product the purchaser accepts the terms and conditions of all applicable Limited Use Label Licenses.

Expression cassette of pASG-IBA44

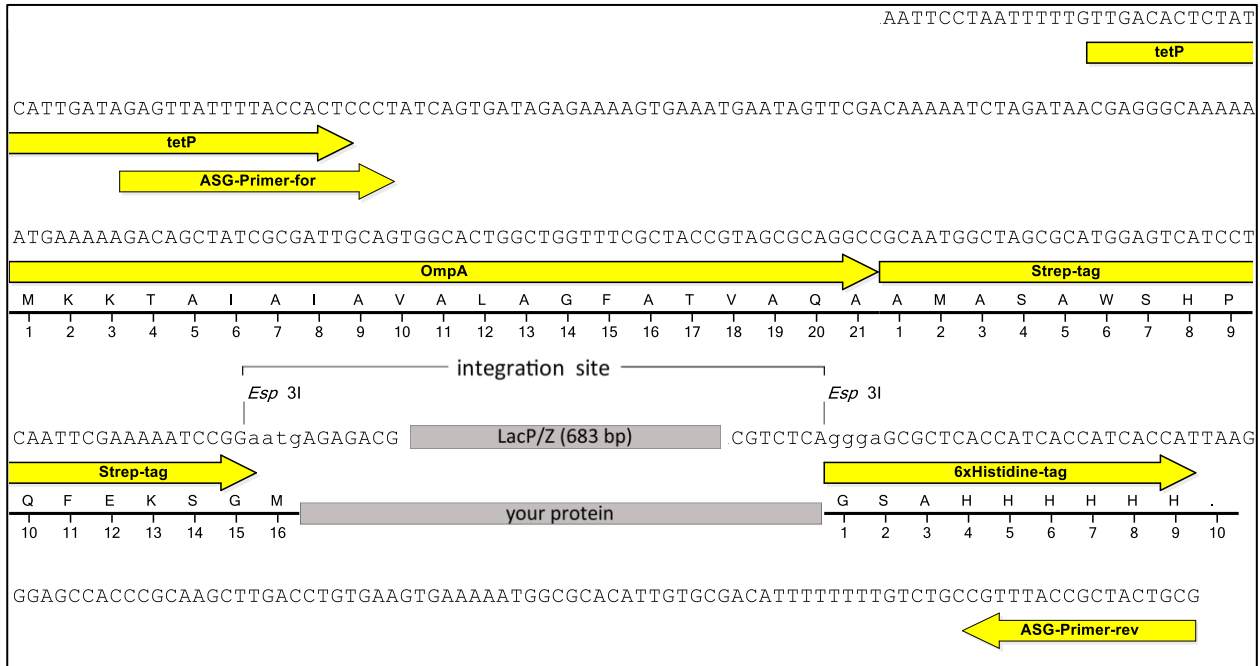

LacP/Z cassette = contains LacZ alpha fragment under control of a separate promoter, which allows alpha complementation of *LacZ* mutations such as *LacZΔM15* as in *E. coli* DH5α or TOP10.

your protein = after StarGate cloning using *Esp3I* your gene of interest will be located here

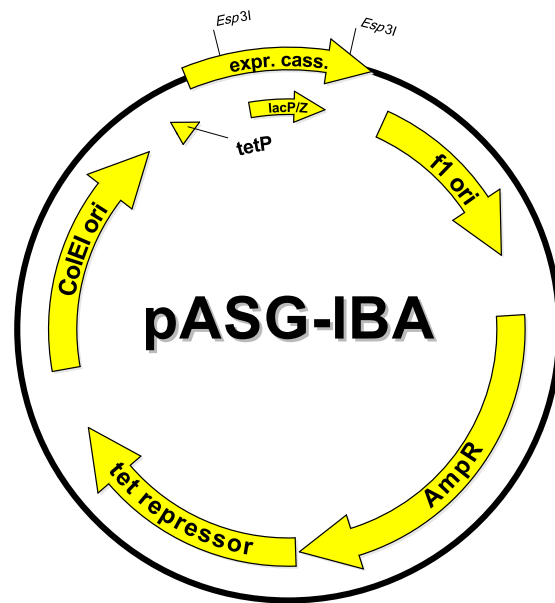

Features	from bp	to bp	Sequencing primer
f1 origin	13	451	ASG-Primer-for 5' - GAGTTATTTTACCACTCCCT -3' ASG-Primer-rev 5' - CGCAGTAGCGGTAAACG -3'
AmpR resistance gene	600	1460	
Tet-repressor	1470	2093	
ColEI ori	2246	2834	
Tet promoter	2939	2975	
forward primer binding site	2959	2978	
OmpA signal sequence	3041	3103	
Strep-tag®II	3104	3148	
LacZ alpha fragment	3377	3778	
6xHistidine-tag	3842	3868	
reverse primer binding site	3942	3958	
total vector length		3958	